


■ MEDICAL SAFETY BOX



PRODUCT PARAMETERS

product name	medical safety box
product color	yellow / white / brown
size	5L / 10L / 20L / other sizes
material	Kraft paper, white board, industrial board, corrugated paper
description	Our products have passed the certification of WHO, and the products can be customized from printing to size

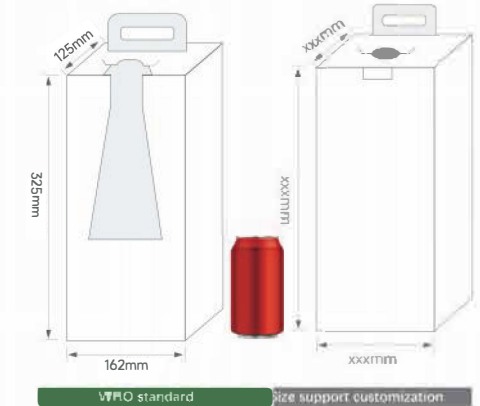

VOLUME: HFD-5L

Colour: white / yellow / brown

size: 325*162*125(mm)/support customization

Capacity : 150 autodisable syringes of 0.5ml with needles



VOLUME: HFD-10L

Colour: white/yellow

size: 325*250*150(mm)/support customization

Capacity : 300 autodisable syringes of 0.5ml with needles



VOLUME: HFD-20L


Colour: white/yellow

size: 400*260*200(mm)/support customization

Capacity : 500 autodisable syringes of 0.5ml with needles

CERTIFIED BY WHO

MATERIAL DISPLAY

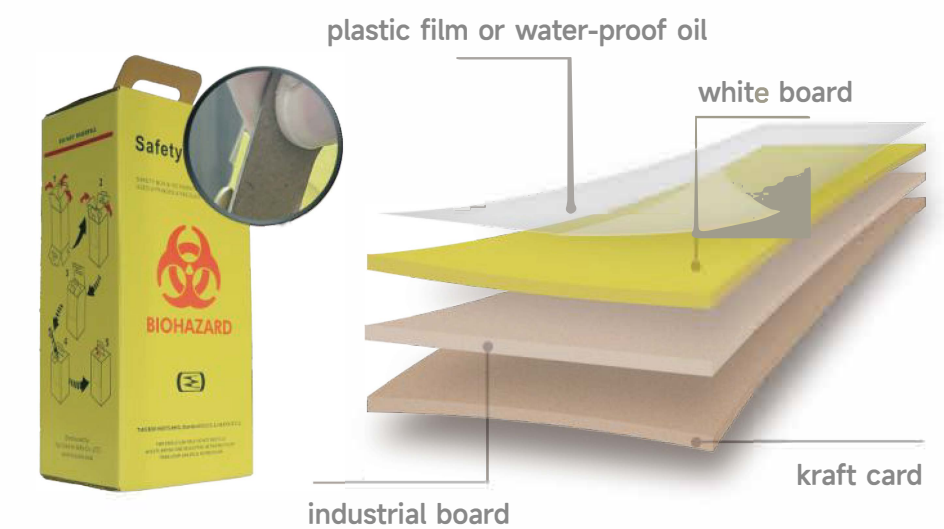

EASY TO ASSEMBLE , EASY TO USE

USAGE METHOD

1. Gently drop the needle and syringe into the opening at the top of the safety box
2. Close the opening of the Safety box after use




Leakproof fluid

Exquisite structure and water-proof coating ensure that the liquid in the box will not leak



Puncture-resistant

Tough composite material ensures no puncture by needles inside the box

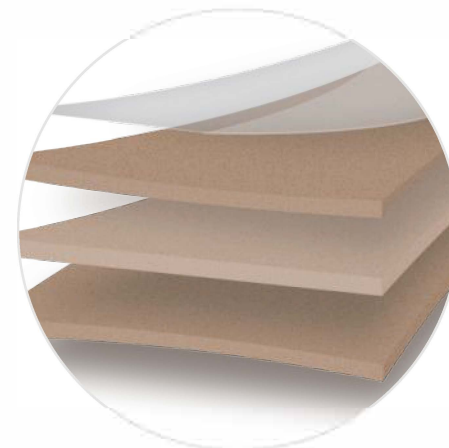

High quality

Specially designed structure even if the product falls from the table, the appearance will not be damaged in any way



High temperature resistance

Even in a high temperature environment of more than 100 degrees, the functionality of the product will not change in any way



Environmentally friendly material

Using environmentally friendly paper materials, the process will not pollute the environment

■ PLASTIC SHARPS CONTAINER

MAKE MOLDS FOR FREE



PRODUCT PARAMETERS


product name	medical plastic sharps container
product color	yellow / red / black / other colours
size	0.5L/1L/1.5L/2L/3L/4.6L/5L/7L/10L/15L/23L ...
material	High-density polyethylene
Instruction	The company can produce sharps containers of any size according to customer requirements and provide free molds





Detachable lid

The diverse design of the hole allows accurate collection of needles, which need to be closed after collection

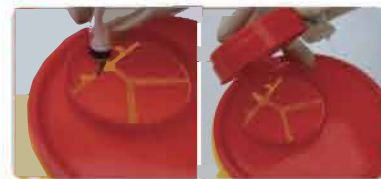


Rotatable lid

The small hole can collect needles, the large hole can collect syringes, various blades, suture needles and so on. When collection is complete, twist the lid to close




Gear lid

The specially designed opening is covered with plastic serrations, which can easily separate the needle from the syringe, and the lid needs to be closed after use

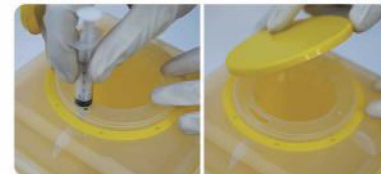


Pushable lid

The specially designed multipurpose opening collects both needles and syringes. Control the lid by pushing and pulling

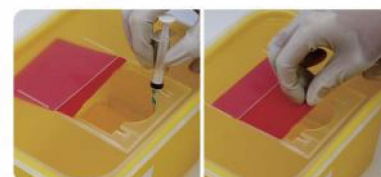


Multifunctional lid-A

The small hole can be used to collect needles or temporarily place syringes. The large hole can collect syringes or other medical supplies. Please close the lid after use

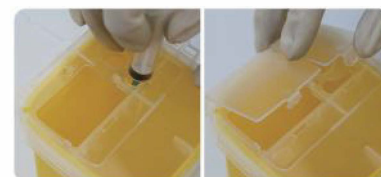

Multifunctional lid-B

Irregular openings accurately collect needles and syringes, and the lid needs to be closed after collection

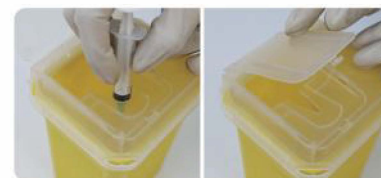


MEDICAL PLASTIC BAGS



PRODUCT PARAMETERS

product name	medical plastic bags
product color	red/yellow/white/ other colours
Size	support customization
material	High-density polyethylene
Features and Application	High portability, good sealing, strong toughness. For medical institutions and non-profit organizations and families.

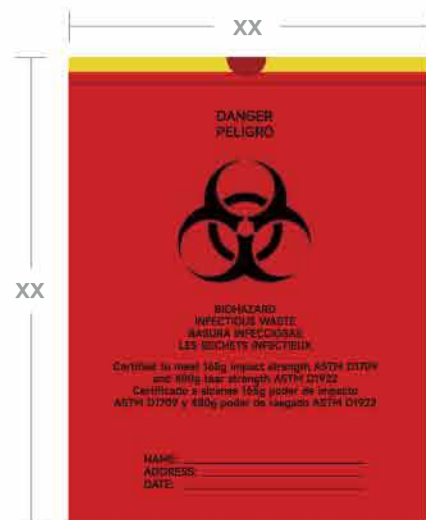

Style : bag with drawstring

color : red yellow and others
print : Support customizatio
size : Support customization
thickness : Support customization



Style : Flat mouth plastic bag

color : red yellow and others
print : Support customizatio
size : Support customization
thickness : Support customization



Style : portable plastic bag

color : red yellow and others
print : Support customizatio
size : Support customization
thickness : Support customization



The company can produce all types of plastic bag products according to customer requirements. Including plastic film, special-shaped bag, medicine bag, etc.

PRODUCT ADVANTAGE



PREVENT LEAKAGE

Plastic bags have high strength and can effectively prevent puncture

Plastic bags are well sealed and liquid will not leak



DEIGN WITH HANDLE

Make it easier to use and seal plastic bags



HIGH TOUGHNESS

Tensile force is not damaged



VARIOUS TYPES

Various types of plastic bags can be customized according to customer requirements.